



APPLIED SCIENTIFIC
INSTRUMENTATION

283 x 110mm Inserts




Eight-Side Insert for MS-9400 with Brass Sliders (I-9408-B)

The I-9408-B slide insert accepts eight 25x75 mm (1" x 3") slides, and has individual spring-loaded brass slider restraints for easy loading and unloading. The unit is recessed to place the bottom of slides about 7.3 mm below the top of the insert. The slides can be easily pre loaded into the insert before installing the insert into the stage. Multiple inserts can be used in this manner to screen a large number of samples.

Depth from top of Insert: 7.3 mm

Overall Thickness: 8.0 mm

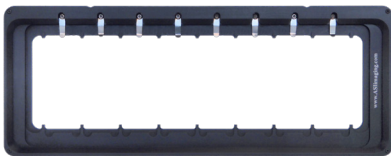

Eight-Side Insert for MS-9400 with Spring Clips (I-9408-C)

The I-9408-C slide insert accepts eight 25x75 mm (1" x 3") slides, and has individual spring clips for secure stability. The unit is recessed to place the bottom of slides about 7.3 mm below the top of the insert. The slides can be easily pre loaded into the insert before installing the insert into the stage. Multiple inserts can be used in this manner to screen a large number of samples.

Depth from top of Insert: 7.3 mm

Overall Thickness: 8.0 mm



60 or 65mm Petri Dish Insert (I-9460, I-9465)

60 or 65mm petri dish insert.

Depth from top of Insert: N/A

Overall Thickness: N/A



2X MULTIWELL (I-9409)

I-9409 2X MV Multiwell slide insert will accept 2 multiwell microplates that are SBS Standard 127.5 x 85 mm (5" x 3 1/3") foot print, for example, the 96-well microplate. The unit places the bottom of a wellplate about 6.7 mm below the top of the insert.



Glass Insert (I-9008)

The I-9008 glass insert provides a large specimen area the full size of the insert. It is made of 4.8 mm (3/16") chemically strengthened soda lime float glass.

Depth from top of Insert: N/A

Overall Thickness: 4.8 mm